

AMP IT UP SPECIAL

USAGE INSTRUCTIONS

MORNING

STEP 1—ANTI-AGE Daily Cleansing Mask

Using wet fingertips, gently and evenly massage a quarter-sized amount onto dry skin. Ensure entire face and neck are covered with a thin coat of product, avoiding contact with eyes. Allow product to dry completely, approximately two minutes. Rinse thoroughly with water; pat dry.

STEP 2—ANTI-AGE Pore Minimizing Toner

Apply evenly to entire face and neck with a gauze pad, avoiding the eye area. Do not rinse off. Allow to dry.

STEP 3—ANTI-AGE Triple Defense Treatment SPF 30

Apply a dime-sized amount over entire face and neck every morning. Avoid eye area.

EVENING

STEP 1—ANTI-AGE Daily Cleansing Mask


Using wet fingertips, gently and evenly massage a quarter-sized amount onto dry skin. Ensure entire face and neck are covered with a thin coat of product, avoiding contact with eyes. Allow product to dry completely, approximately two minutes. Rinse thoroughly with water; pat dry.

STEP 2—ANTI-AGE Pore Minimizing Toner

Apply evenly to entire face and neck with a gauze pad, avoiding the eye area. Do not rinse off. Allow to dry.

STEP 3—AMP MD™ Micro-Channeling Tool

Use at least three nights per week, increasing to nightly use as tolerated. With moderate and comfortable pressure, roll AMP MD over entire face and neck 4-10 times. For best results, divide the face into sections (forehead, right side, left side, chin, neck) and treat each section before continuing onto the next. Change the orientation of rolling over each section (as illustrated) to ensure complete surface coverage.

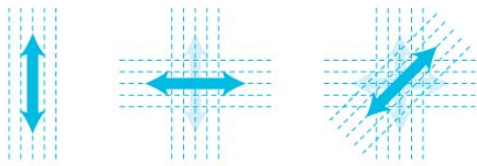

After each use, rinse AMP MD with very warm running water. Twice per week, clean thoroughly using AMP MD Purification Tablets, using a new tablet for each use:

Fill Cleansing Vial with very warm (not hot) water to marked fill line and drop one Purification Tablet into water.

Securely snap the AMP MD handle into the cap of the cleansing vial and immediately immerse the AMP MD into the cleansing solution.

When blue color fades, approximately 15 minutes, remove AMP MD and rinse thoroughly under very warm running water. Do not allow AMP MD to remain in purification solution overnight.

Respect the Micro-Needles: AMP MD is constructed of the highest quality surgical-grade acupuncture needles for maximum efficacy. Always store AMP MD in the clear storage case or snapped into the cap of the Cleansing Vial, and do not place AMP MD roller head directly on hard surfaces or counter tops.

For your safety: AMP MD should not be used on the eyes, lips, irritated or infected skin, and is not recommended for use on hyper-sensitive skin. Keep out of reach of children.

STEP 4—ANTI-AGE Night Renewing Serum

Open one capsule and apply serum evenly to entire face and neck, avoiding the eye area.

STEP 5—ANTI-AGE Overnight Restorative Cream

Apply a dime-sized amount evenly over the entire face and neck every evening, avoiding the eye area.